

## MONTGOMERY® LARYNGEAL KEEL

### DESCRIPTION:

The **Montgomery® Laryngeal Keel**, with its unique umbrella shape, has proven useful in the repair of anterior subglottic stenosis (web), with or without vocal cord paralysis. Use of the keel can be made following repair of stenosis, or; if use of a Montgomery® Laryngeal Stent is indicated, following removal of the stent to aid in formation of a sharp anterior commissure.

The keel also has application following hemilaryngectomy to prevent stenosis. The intralaryngeal extension of the keel is thinner than the umbrella of the extralaryngeal portion. The extralaryngeal portion, when secured to thyroid laminae with a figure-of-eight suture, covers the thyrotomy repair with its umbrella. Tightly secured, the keel prevents re-stenosis or the formation of granulation tissue.

The keel is made of soft medical grade silicone (clear or radiopaque) and is available in three sizes.

### INDICATIONS FOR USE:

- Repair of anterior glottic stenosis.
- Formation of a sharp anterior commissure following laryngeal stent removal.
- Prevention of stenosis following hemilaryngectomy.



### SUPPLIED:

Sterile, 1 keel per carton

### PRODUCT CODES:

<b>323040</b>	<b>SIZE 12mm KEEL, CLEAR</b>
<b>323040R</b>	<b>SIZE 12mm KEEL, RADIOPAQUE</b>
<b>323050</b>	<b>SIZE 14mm KEEL, CLEAR</b>
<b>323050R</b>	<b>SIZE 14mm KEEL, RADIOPAQUE</b>
<b>323060</b>	<b>SIZE 16mm KEEL, CLEAR</b>
<b>323060R</b>	<b>SIZE 16mm KEEL, RADIOPAQUE</b>

Dimensions (mm)	A
Size 12 mm	12
Size 14 mm	14
Size 16 mm	16

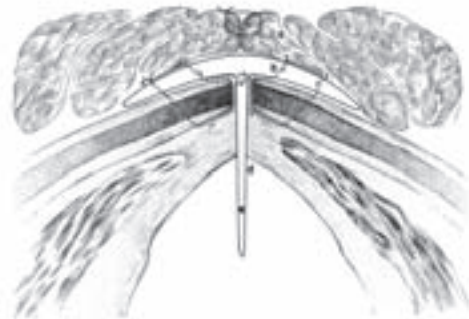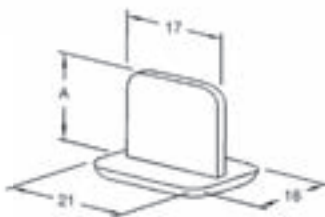

### FEATURES

- Implant grade silicone.
- Non-adherent smooth surface.
- Available in clear or radiopaque material.

